


# Glycine-rich peptides from fermented *Chenopodium formosanum* sprout as an antioxidant to modulate the oxidative stress

Chen-Che Hsieh , Shu-Han Yu , Hsing-Chun Kuo , Darin Khumsupan , Hsiao-Chu Huang , Yu-Wei Liou , Chen-Yu Kao , Szu-Chuan Shen , Kuan-Chen Cheng

## Abstract

*Rhizopus oligosporus* was used to ferment *Chenopodium formosanum* sprouts (FCS), yielding a glycine-rich peptide (GGGGGKP) from <2 kDa proteins. This peptide showed antioxidant activity and was confirmed to be safe, non-toxic, and non-allergenic via SwissADME and other tools. In silico and in vitro studies revealed that it inhibits v-rel phosphorylation in the NF-κB pathway, reducing oxidative stress and inflammation, indicating potential biomedical applications.

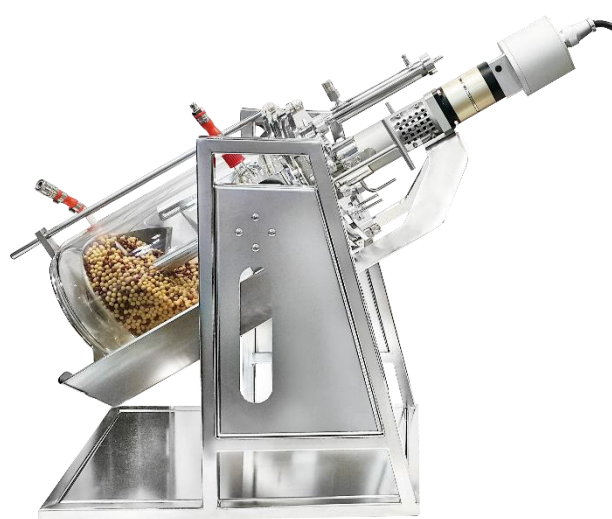

**Winpact Model: FS-V-SA05P**

## Introduction

Fermentation is known to enhance the bioactivity of plant-based proteins. *Chenopodium formosanum* sprouts are nutrient-rich and contain bioactive compounds. Prior research has highlighted their antioxidant and anti-inflammatory potential. This study focuses on isolating <2 kDa peptides, identifying glycine-rich sequences (GGGGGKP), and evaluating their biological properties and safety through both computational and experimental assays.

### US Office

19959 Sea Gull Way  
Saratoga, CA 95070  
U.S.A.  
T/ +1-408-366-9866  
F/ +1-408-446-1107

### Taiwan Office

No. 156, Sec. 1, Guoji Rd., Taoyuan Dist.,  
Taoyuan City 330041, Taiwan  
T/+886-3-3762878  
F/+886-3-3761310

### Shanghai Office

Room 2519, No. 306 Taopu Road,  
Putuo District, Shanghai, China  
Postal code: 200333  
Cell:+86-13621808204  
T/ +86-21-50795277  
F/ +86-21-50795277

### India Office

D.No.7-143, 2nd Floor,  
St.No.2,Nagendra Nagar,  
Habsiguda, Hyderabad-500007,  
India  
T/ +91-40-27001515  
T/ +91-40-27001586

## Materials and Methods

### *Fermentation Process:*

1.5 kg of *Chenopodium formosanum* sprouts were inoculated with 75 mL of *Rhizopus oligosporus* spore suspension ( $10^6$  spores/mL). Solid-state fermentation was carried out in a \*\*bioreactor model FSeV-SA05P (Winpack FS-V-SA05P)\*\* at 35°C, with aeration at 0.4 vvm and rotation at 5 rpm, for 4 days.

## Results

- Bioreactor fermentation (BF) produced significantly higher free peptide content and enzyme activities compared to plate fermentation (PF).
- Identified peptides potentially act as DPP-IV and ACE inhibitors.
- Over twenty new metabolites, including aromatics, amines, fatty acids, and carboxylic acids, were generated specifically in BF.
- Enhanced antioxidant capacity was demonstrated in BF-derived products..

## References

Glycine-rich peptides from fermented *Chenopodium formosanum* sprout as an antioxidant to modulate the oxidative stress

<https://doi.org/10.1016/j.foodres.2023.112707>

### US Office

19959 Sea Gull Way  
Saratoga, CA 95070  
U.S.A.  
T/ +1-408-366-9866  
F/ +1-408-446-1107

### Taiwan Office

No. 156, Sec. 1, Guoji Rd., Taoyuan Dist.,  
Taoyuan City 330041, Taiwan  
T/+886-3-3762878  
F/+886-3-3761310

### Shanghai Office

Room 2519, No. 306 Taopu Road,  
Putuo District, Shanghai, China  
Postal code: 200333  
Cell: +86-13621808204  
T/ +86-21-50795277  
F/ +86-21-50795277

### India Office

D.No.7-143, 2nd Floor,  
St.No.2,Nagendra Nagar,  
Habsiguda, Hyderabad-500007,  
India  
T/ +91-40-27001515  
T/ +91-40-27001586

major  
science  
[www.majorsci.com](http://www.majorsci.com)  
[info@majorsci.com](mailto:info@majorsci.com)