

Winpact | Winpact Shaking Incubator



CE SI-200



Polycarbonate Erlenmeyer Flask
(Baffled Base)


Suitable with SI series products



Baffled Base



Individually Wrapped

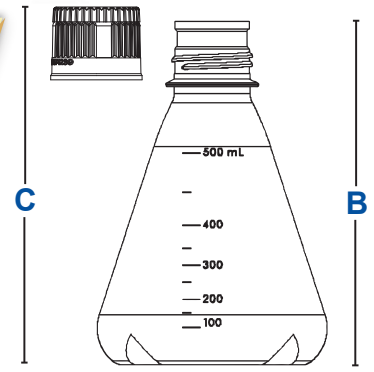


A • Base Diameter
B • Height (w/o Cap)
C • Height (w/ Cap)



Patented DuoCAP® System



Standard with all TriForest Shaker Flasks, it allows aerobic/anaerobic culture of microbes and algae without replacing the cap. The 0.22 um PTEE filter allows sterile exchange of gases and air.

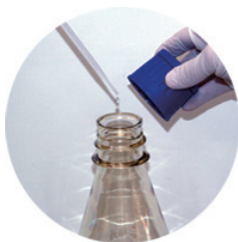

Add Culture Media & Inoculum



Add DuoCAP™ for Anaerobic Culture



Remove dustguard cover for Aerobic Culture



Move to Shaker Incubator

Specification

Baffled Base Cat.	FBC2000S	FBC1000S	FBC0500S	FBC0250S	FBC0125S
Size / Volume	2000 mL	1000 mL	500 mL	250 mL	125 mL
Closure	53B	53B	38/430	38/430	38/430
Cap Type	DuoCAP®	DuoCAP®	DuoCAP®	DuoCAP®	DuoCAP®
Units / Pack	6	6	12	12	24
Packs/Case (Total Units)	4(24)	4(24)	4(48)	6(72)	6(144)

